

ARJOHUNTLEIGH

GETINGE GROUP

MINUET 2

ADVANCED LONG TERM CARE BED



...with people in mind

In a comparative study, the *Minuet*® 2 with *Pro-Contour*® Advance “gave statistically significant lower interface pressure profiles and exerted lower pressure gradients between the flat and Fowler’s position”¹

“The design of the mattress platform, until recently of lesser interest, can have a significant effect on pressure exerted at the mattress body interface”¹

Intelligent Profiling System

To prevent the potential for the resident to slide down the bed, upon lowering of the backrest, the knee-break will stay raised until the backrest is at the same angle (15 degrees) as the knee-break. At that time, both the backrest and knee-break will descend together until both reach a flat position simultaneously which allows the resident greater stability during repositioning of the *Minuet*® 2 frame.


This is unique - on other beds the knee-break descends regardless of the backrest angle, which can result in the resident sliding down the bed.

The intelligent profiling system maximizes comfort, improves postural support and helps to reduce pressure and shear which can help to minimize the risk of pressure ulcers. It also can minimize moving and handling risks for the caregiver as the resident is less prone to slide down the bed.


Pro-Contour® Advance Profiling

Pro-Contour Advance Profiling means that the *Minuet* 2 is anatomically designed to expand the resident’s sacral area on the frame as the head of the bed is raised. This provides reduced pressure in the sacral area and increased comfort for the resident. This offers major benefits over conventional profiling beds:

- Enhanced pressure and shear reduction
- Excellent postural support
- Exceptional comfort

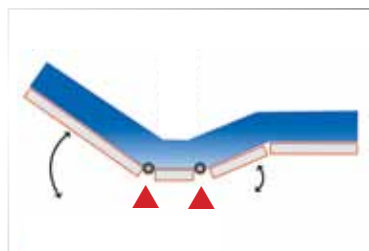

▲ Conventional bed pivot points

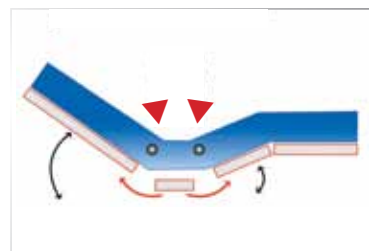

▲ *Pro-Contour*® articulation mimics the resident's natural back and thigh pivot points



Optional Siderail



Half length mobility assist siderail

CM-ACC23-1

Auto-Contour Button



On the *Minuet 2* handset, the orange Auto-Contour button allows:

- *Pro-Contour Advance* profiling
- Effortless upright positioning
- Intelligent profiling system

All of the handset buttons are raised and embossed for ease of use

Lifting Poles



Lifting Pole accessory

CM-ACC09-1
Single position

CM-ACC10-1
Multi position

Mattress Extension



An 5" (13 cm) mattress extension is available

Pump Bracket



Pump Bracket accessory

CM-ACC15-1

Product Specifications

Minuet 2 Bed	
Length	88½" (225 cm)
Width (without siderails)	37" (94 cm)
Width (with siderails)	41¼" (105 cm)
In-bed length	80" (203 cm)
Height range	8 ¾"-25 ¼" (25-66 cm)
Backrest angle (Maximum)	73°
Knee-break elevation (Maximum)	32°
Safe Working Load (SWL)	396 lbs (180 kg)
Weight	195.7 lbs (88.9 kg)



Key Advantages

- Battery back-up allows articulation of the bed in the event of a power failure
- *Pro-Contour* Advance Profiling enhances the reduction of pressure, shearing and friction
- Wide range of active and reactive pressure redistributing support surfaces

Features	
Electric height & backrest adjustment	standard
Electric knee-break adjustment	standard
Battery back-up	standard
Handset with raised, embossed buttons	standard
Flexible handset holder (CM-ACC13-1)	optional
Raising, retracting backrest	standard
<i>Pro-Contour</i> Advance Profiling	standard
Auto-Contour button	standard
Ratchet calf section adjustment	standard
Integral bed extension with mattress retainer	standard
Beech effect head and foot end panels	standard
Slatted mattress platform	standard
Kick stop brakes	standard
Disassembles into 4 sections for transport/storage	standard

¹ Call E, Baker B. How does bed frame design influence tissue interface pressure? A comparison of four different technologies designed for long-term or home care. Journal of Tissue Viability. 2007; 17:22-29

Use only ArjoHuntleigh designed parts, i.e. parts specifically designed for the purpose, on equipment and products supplied by ArjoHuntleigh. As our policy is one of continuous development, we reserve the right to modify designs and specifications without prior notice. ® and ™ are trademarks belonging to the ArjoHuntleigh group of companies. © ArjoHuntleigh, 2010

Printed in USA • minuet2_broch_033110

ARJOHUNTLEIGH

GETINGE GROUP

www.ArjoHuntleigh.com

GETINGE GROUP is a leading global provider of products and systems that contribute to quality enhancement and cost efficiency within healthcare and life sciences. We operate under the three brands of ArjoHuntleigh, GETINGE and MAQUET. **ArjoHuntleigh** focuses on patient mobility and wound management solutions. **GETINGE** provides solutions for infection control within healthcare and contamination prevention within life sciences. **MAQUET** specializes in solutions, therapies and products for surgical interventions and intensive care.

ArjoHuntleigh Inc.
2349 West Lake Street
Addison, Illinois 60101
Phone: 800-323-1245
Fax: 630-925-7969

ArjoHuntleigh Canada Inc.
1575 South Gateway Road Unit C
Mississauga, Ontario L4W 5J1
CANADA
Phone: 800-665-4831
Info.Canada@ArjoHuntleigh.com

